-  (EN) Soy Protein Rapid Kit
-  (FR) Kit Rapide Protéine Soja
-  (DE) Soja Protein Rapid Kit
-  (IT) Kit rapido per la rilevazione delle proteine di soia
-  (ES) Kit Rápido para Proteína de Soya
-  (NL) Soja Proteïne Sneltest
-  (SV) Soy Protein Rapid Kit
-  (DA) Soya Protein Rapid Kit
-  (NO) Soyaprotein hurtigkit
-  (PT) Kit Rápido para Proteína de Soja
-  (EL) Κιτ Ταχείας Δοκιμασίας Πρωτεΐνης Σόγιας
-  (PL) Szybki test na wykrywanie białek obecnych w soi
-  (RU) Rapid набор (протеин сои)
-  (TR) Soya Proteinini Hızlı Test Kiti
-  (JA) 大豆プロテイン迅速測定キット
-  (ZH) 大豆蛋白快速检测试剂盒
-  (TH) ชุดทดสอบโปรตีนของถั่วเหลืองแบบรวดเร็ว
-  (KO) 콩 단백질 신속 키트
-  (ID) Rapid Kit Protein Kacang Kedelai
-  (AR) طقم سريع بروتين فول الصويا

# Product Instructions

## Soy Protein Rapid Kit

Lateral Flow Devices (LFD) for qualitative analysis of soy proteins.

### Product Description and Intended Use

The Neogen® Soy Protein Rapid Kit is intended for screening for the presence of soy proteins in clean-in-place (CIP) final rinse water, environmental swab samples, food ingredients and processed food products.

The Neogen Soy Protein Rapid Kit utilizes a lateral flow device (LFD) that is an immunochromatographic test method utilizing antibodies which are specific for the detection of soy proteins. Positive results are visualized by the presence of three lines: a control line, a hook line and a test line, when the soy protein is present in the range of concentrations from a lower level of 2 ppm and a higher level of approximately 2 mg/mL of soy protein in saline solution (0.9% sodium chloride). These limits may vary depending on the matrix tested. Samples containing more than 1% soy may result in an invalid test which will result in the absence of the hook line and/or test line.

The Neogen Soy Protein Rapid Kit is intended for use in the food and beverage industry by trained personnel. Neogen has not documented the use of this product in industries other than food or beverage. For example, Neogen has not documented this product for testing pharmaceutical, cosmetic, clinical or veterinary samples. The Neogen Soy Protein Rapid Kit has not been evaluated with all possible food products, food processes and testing protocols.

The Neogen Soy Protein Rapid Kit contains 25 tests, described in Table 1.

**Table 1.** Kit components

Item	Identification	Quantity	Storage
Neogen® Soy Protein Lateral Flow Device (LFD)	Lateral flow device in a plastic cassette	25 devices individually packed.	Store at 2-8°C. Do not freeze.
Neogen® Extraction Buffer	Bottle with Extraction Buffer	1 bottle with 50 mL	Store at 2-8°C. Do not freeze.
Dilution Tubes	Microcentrifuge tube (2.2 mL volume capacity)	26 tubes	Store in a clean dry place.

Materials not provided in the kit:

- Swabs and pipettes.
- The use of vortex, timer and balance are recommended but not required for all samples.
- The use of a centrifuge is required for all chocolate and gum samples and it is recommended, but not required for all solid samples.

### Safety

The user should read, understand, and follow all safety information in the instructions for the Neogen Soy Protein Rapid Kit. Retain the safety instructions for future reference.

**⚠ WARNING:** Indicates a hazardous situation, which, if not avoided, could result in death or serious injury and/or property damage.

**NOTICE:** Indicates a potentially hazardous situation, which, if not avoided, could result in property damage.

### **⚠ WARNING**

**To reduce the risks associated with inaccurate results:**

- Neogen has not documented the use of Neogen Soy Protein Rapid Kit in industries other than food or beverage. For example, Neogen has not documented this product for testing pharmaceutical, cosmetic, clinical or veterinary samples.
- Neogen Soy Protein LFD should be read  $11 \pm 1$  minutes after sample has been loaded on the lateral flow device.
- The Neogen Extraction Buffer is designated for use with a specific lot of Neogen Soy Protein LFD. Do NOT interchange Neogen Soy Protein Rapid Kit components with other lots or kits.
- The Neogen Extraction Buffer is designated for use with a specific lot of Neogen Soy Protein LFD. Dispose of any remaining Neogen Extraction Buffer once all Neogen Soy Protein Lateral Flow Devices have been used.

- Store the Neogen Soy Protein Rapid Kit as indicated on the package and in the product instructions.
- Always use the Neogen Soy Protein Rapid Kit by the expiration date.
- Always use the Neogen Soy Protein Rapid Kit at 20-25°C temperature.
- Neogen® Allergen Protein Testing Kits are not intended for the detection of hydrolyzed proteins.
- The Neogen Allergen Protein Testing Kits are designed to detect proteins from processed food once they are solubilized into Neogen Extraction Buffer. Some food processing methods may limit detection of these target proteins.
- Some food processing may affect the detection of food proteins with Neogen Allergen Protein Testing Kits. Users should verify that the method is fit for purpose to meet user's requirements.

**To reduce the risks associated to false negative results:**

- Use the Neogen Soy Protein Rapid Kit for food and environmental samples that have been validated internally or by a third party.

**To reduce the risks associated with exposure to chemicals:**

- Neogen Soy Protein Rapid Kit is intended for use in the food and beverage industries by trained personnel.

**NOTICE****To reduce the risks of inaccurate results:**

- Refer to the Interpretation of Results section of the product instructions, to ensure accurate interpretation of the Neogen Soy Protein LFD.

Consult the Safety Data Sheet for additional information.

For information on documentation of product performance, visit our website at [www.neogen.com](http://www.neogen.com) or contact your local Neogen representative or distributor.

**User Responsibility**

Users are responsible for familiarizing themselves with product instructions and information. Visit our website at [www.neogen.com](http://www.neogen.com), or contact your local Neogen representative or distributor for more information.

**As with all test methods used for food analysis the test matrix can influence the results.** When selecting a test method, it is important to recognize that external factors such as sampling methods, testing protocols, sample preparation, handling, and laboratory technique may influence results. The food sample itself may influence results.

It is the user's responsibility in selecting any test method or product to evaluate a sufficient number of samples to satisfy the user that the chosen test method meets the user's criteria.

It is also the user's responsibility to determine that any test methods and results meet its customers' and suppliers' requirements.

As with any test method, results obtained from use of any Neogen Food Safety product do not constitute a guarantee of the quality of the matrices or processes tested.

**Limitation of Warranties/Limited Remedy**

EXCEPT AS EXPRESSLY STATED IN A LIMITED WARRANTY SECTION OF INDIVIDUAL PRODUCT PACKAGING, NEOGEN DISCLAIMS ALL EXPRESS AND IMPLIED WARRANTIES, INCLUDING BUT NOT LIMITED TO, ANY WARRANTIES OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR USE. If any Neogen Food Safety Product is defective, Neogen or its authorized distributor will, at its option, replace or refund the purchase price of the product. These are your exclusive remedies. Please contact your Neogen representative or authorized Neogen distributor for any further questions.

**Limitation of Neogen Liability**

NEOGEN WILL NOT BE LIABLE FOR ANY LOSS OR DAMAGES, WHETHER DIRECT, INDIRECT, SPECIAL, INCIDENTAL OR CONSEQUENTIAL DAMAGES, INCLUDING BUT NOT LIMITED TO LOST PROFITS. In no event shall Neogen's liability under any legal theory exceed the purchase price of the product alleged to be defective.

**Storage and Disposal**

Store all Neogen Soy Protein Rapid Kit components at 2-8°C.

Neogen Soy Protein Rapid Kit components should not be frozen, exposed to UV light or exposed to prolonged heat (>30°C).

Neogen Soy Protein Rapid Kit components should not be used past the expiration date. Expiration date and lot number are noted on the outside label of the box.

Please note that each Neogen Extraction Buffer lot is validated specifically for each LFD lot and is not interchangeable with any other lots or kits.

Dispose according to current local/regional/industry standards and regulations.

## Instructions for Use

Follow all instructions carefully. Failure to do so may lead to inaccurate results. Ensure that all Neogen Soy Protein Rapid Kit components are at ambient temperature (20-25°C) prior to use.


## Sample Analysis

### 1. CIP Final Rinse Water Samples

- 1.1 Label one microcentrifuge tube for each CIP sample.
- 1.2 Add 800  $\mu$ L of Neogen Extraction Buffer to a labeled microcentrifuge tube.
- 1.3 Add 200  $\mu$ L of CIP final rinse water sample. Shake vigorously or vortex for 15 seconds to mix thoroughly to obtain an extracted sample.


Note: The pH of extracted sample should be between 5 and 10. Proceed to Troubleshooting section for additional information.

- 1.4 Remove one Neogen Soy Protein LFD from package and place on a clean, dry, flat surface.
- 1.5 Transfer 100  $\mu$ L of the extracted sample prepared in 1.3 using a clean pipette or pipette tip and apply it to the sample well on the Neogen Soy Protein LFD. Start the timer for 11  $\pm$  1 minutes. Proceed to Interpretation of Results section.
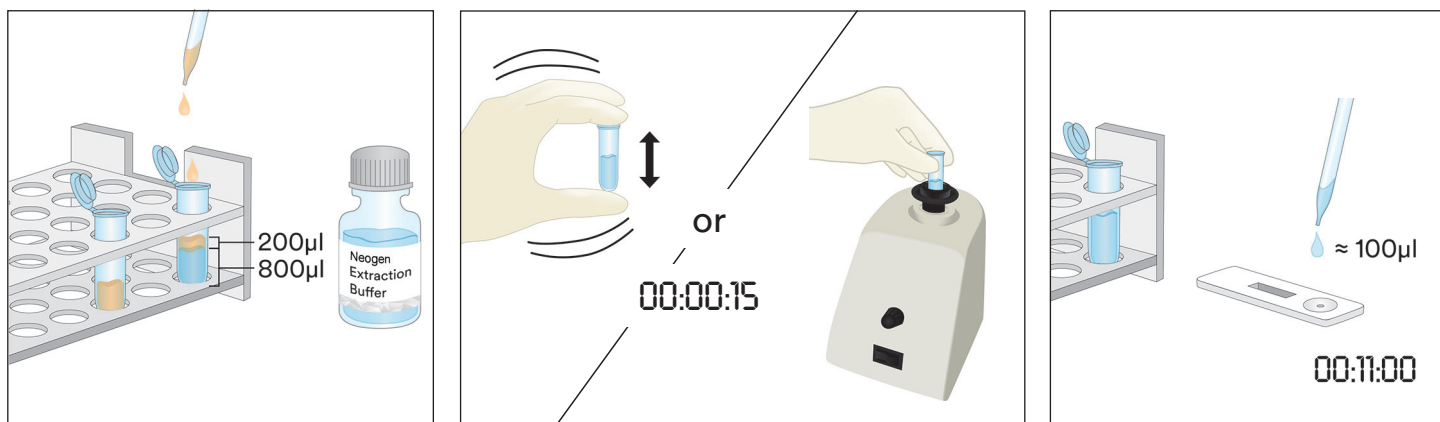
### 2. Environmental Swab Samples

- 2.1 Label one microcentrifuge tube for each Environmental Swab Sample.
- 2.2 Add 500  $\mu$ L of Neogen Extraction Buffer into a labeled microcentrifuge tube.
- 2.3 Take a clean swab and dip the entire tip into the microcentrifuge tube wetting the tip with Neogen Extraction Buffer. Gently express excess liquid from tip by pressing swab tip lightly on the inside of the tube.
- 2.4 Take wetted swab and survey a 10 X 10 cm surface area maintaining the swab at a 30° angle with the surface. Rub the swab slowly and thoroughly over the surface area. Rub the swab three times over this surface, reversing direction between alternating strokes.
- 2.5 Take swab and insert it back into the pre-labeled tube and swirl swab several times to release any residues that might be on the surface of swab into the Neogen Extraction Buffer. Break off swab tip in tube, cap tightly and mix well to obtain an extracted sample.
- 2.6 Remove one Neogen Soy Protein LFD from the package and place on a clean, dry, flat surface.
- 2.7 Transfer 100  $\mu$ L of the extracted sample prepared in 2.5 using a clean pipette or pipette tip and apply it to the sample well on the Neogen Soy Protein LFD. Start the timer for 11  $\pm$  1 minutes. Proceed to Interpretation of Results section.

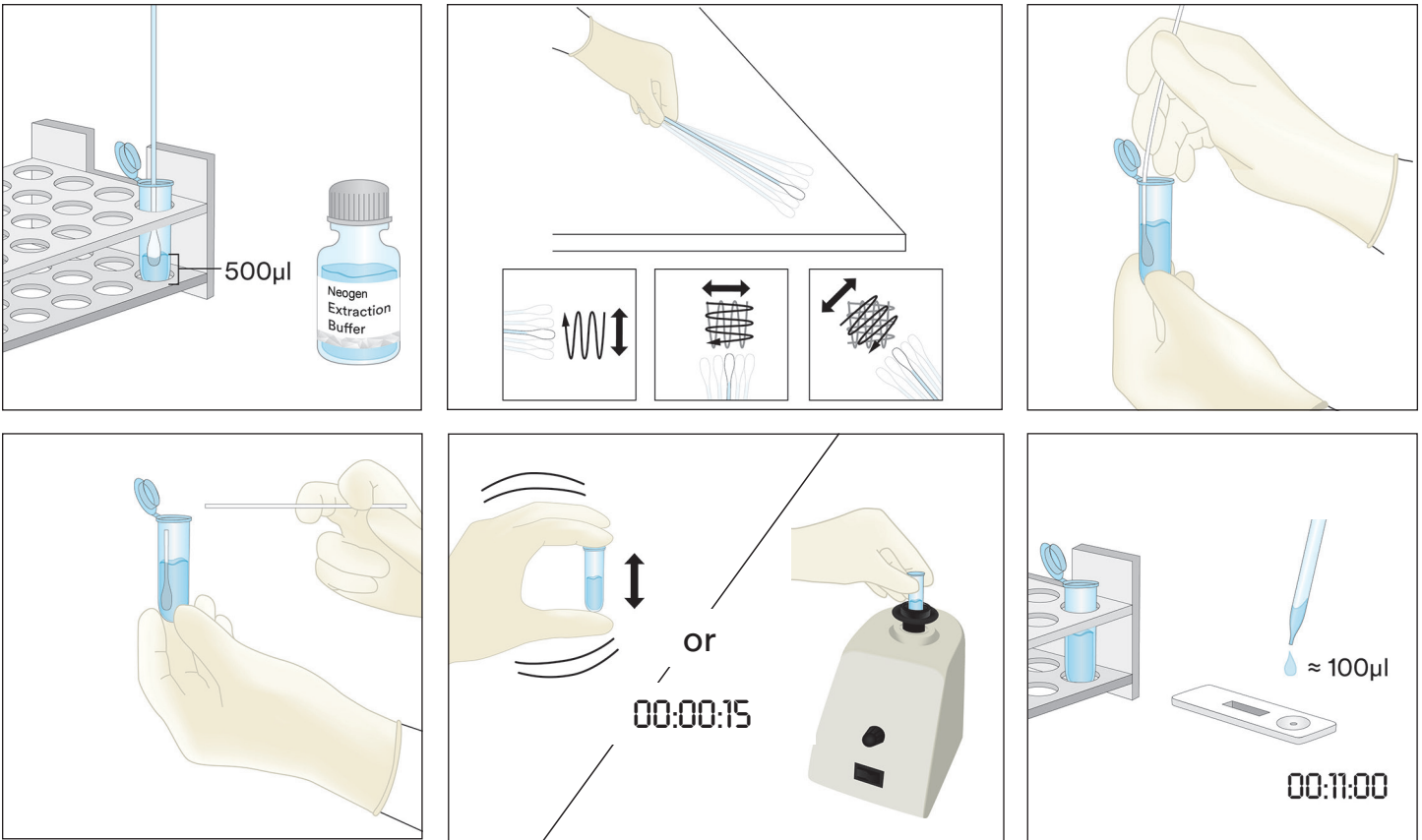
### 3. Liquid Samples Except Liquid Chocolate Samples

- 3.1 Label one microcentrifuge tube for each Liquid sample.
- 3.2 Measure 900 µL of Neogen Extraction Buffer into a labeled microcentrifuge tube.
- 3.3 Add 100 µL of a well-mixed sample. Shake vigorously or vortex for 15 seconds to mix thoroughly to obtain an extracted sample.


Note: The pH of extracted sample should be between 5 and 10. Proceed to Troubleshooting section for additional information.

- 3.4 Remove one Neogen Soy Protein LFD from package and place on a clean, dry, flat surface.
- 3.5 Transfer 100 µL of the extracted sample prepared in 3.3 from the middle (aqueous) layer using a clean pipette or pipette tip and apply it to the sample well on the Neogen Soy Protein LFD. Start the timer for  $11 \pm 1$  minutes. Proceed to Interpretation of Results section.
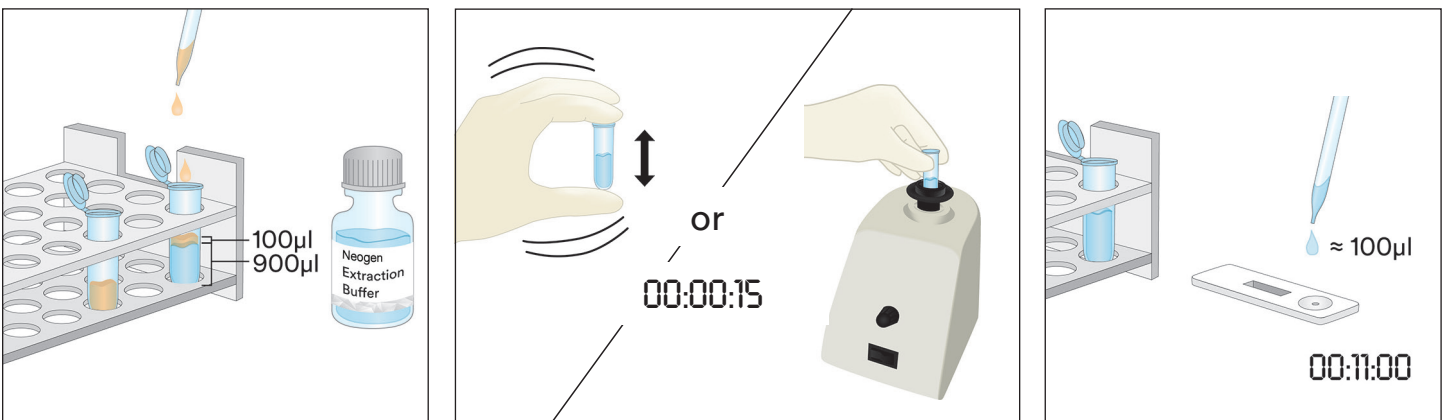
#### 4. Liquid Chocolate Samples

- 4.1 Label two microcentrifuge tubes for each Liquid Chocolate Sample.
- 4.2 To make a prepared sample, add 500  $\mu\text{L}$  of a well-mixed Liquid Chocolate Sample and add 500  $\mu\text{L}$  of Neogen Extraction Buffer into one labeled microcentrifuge tube and shake vigorously to mix thoroughly or vortex for approximately 15 seconds.
- 4.3 Add 900  $\mu\text{L}$  of Neogen Extraction Buffer into the second labeled microcentrifuge tube and add 100  $\mu\text{L}$  of the Prepared Sample from step 4.2. Shake vigorously to mix thoroughly or vortex for approximately 15 seconds.
- 4.4 Centrifuge for 20-30 seconds at 5000-7000 rpm (3000 x g). The supernatant is the extracted sample.
- 4.5 Remove one Neogen Soy Protein LFD from package and place on a clean, dry, flat surface.
- 4.6 Transfer 100  $\mu\text{L}$  of the extracted sample prepared in 4.4 from the middle (aqueous) layer using a clean pipette or pipette tip and apply it to the sample well on the Neogen Soy Protein LFD. Start the timer for  $11 \pm 1$  minutes. Proceed to Interpretation of Results section.
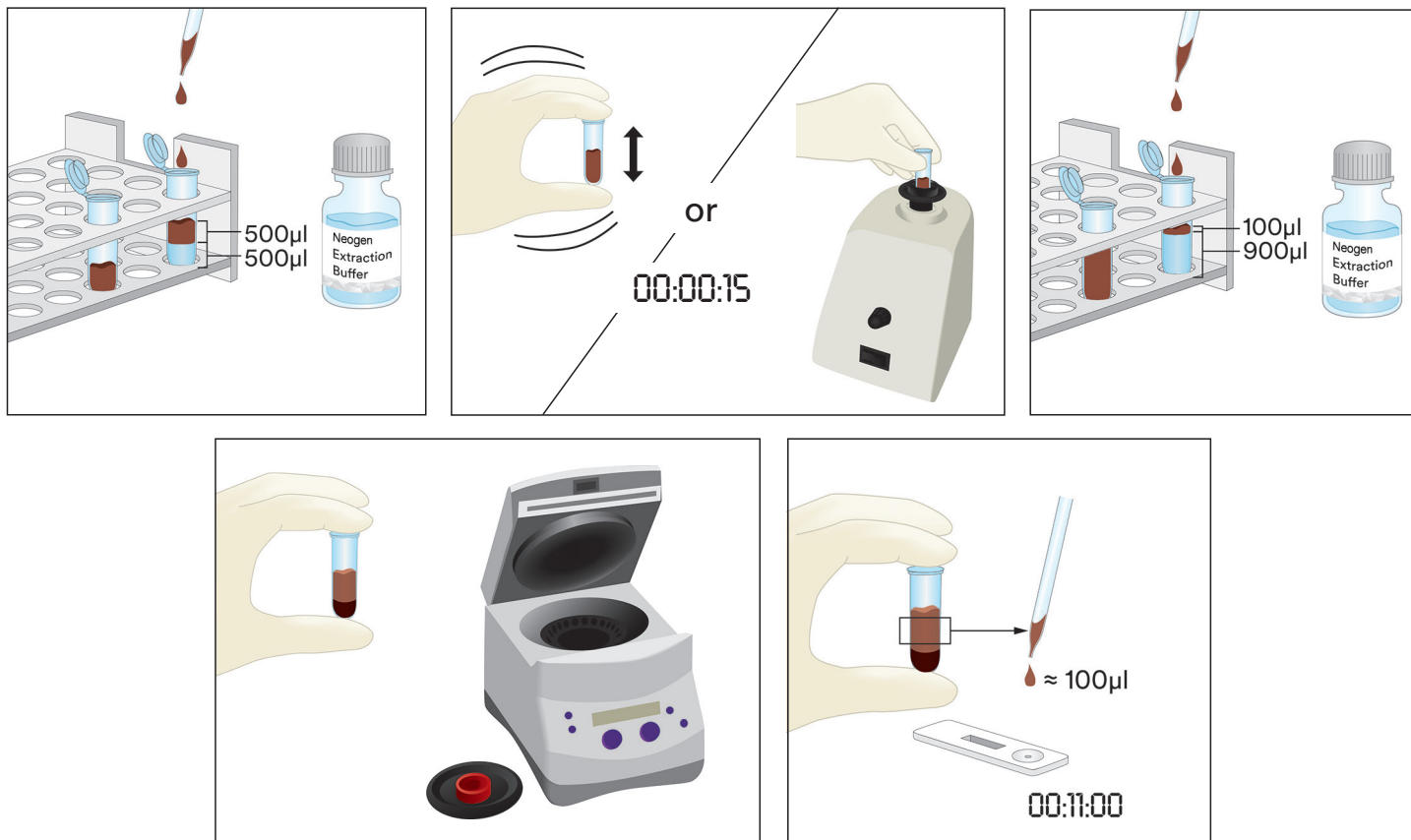
#### 5. Solid Samples

- 5.1 Label one microcentrifuge tube for each non-liquid sample.
- 5.2 Grind sample into a fine homogeneous powder.
- 5.3 Measure 0.2 g of sample into a labeled microcentrifuge tube.
- 5.4 Add 1.8 mL of the Neogen Extraction Buffer to the sample in the microcentrifuge tube. Shake vigorously or vortex for 15 seconds to mix thoroughly to obtain an extracted sample.  
Note: Melt chocolate samples.
- 5.5 Centrifuge for 20-30 seconds at 5000-7000 rpm (3000 x g). The supernatant is the extracted sample.
- 5.6 Remove one Neogen Soy Protein LFD from package and place on a clean, dry, flat surface.
- 5.7 Transfer 100  $\mu\text{L}$  of the extracted sample prepared in 5.5 from the middle (aqueous) layer using a clean pipette or pipette tip and apply it to the sample well on the Neogen Soy Protein LFD. Start the timer for  $11 \pm 1$  minutes. Proceed to Interpretation of Results section.


### Interpretation of Results

The control line is next to the letter C on the Neogen Soy Protein LFD. The test line is next to the letter T on the Neogen Soy Protein LFD. The hook line is next to the letter H and/or in between the T and C lines. Note: Some LFD cassettes might not have an H marking on the side of the plastic cassette.

Read LFD at 11 ± 1 minutes after application of the sample. A sample is considered to be:

- Negative for soy protein when only the two lines furthest from the well, the control and the hook line, are visible on the Neogen Soy Protein LFD.
- Positive for soy protein when all three lines (test, hook and control line) are visible on Neogen Soy Protein LFD.
- Invalid, if the Neogen Soy Protein LFD does not develop the middle hook line. The sample may be overloaded with protein and needs to be diluted.
- Invalid, if the Neogen Soy Protein LFD does not develop the control line.

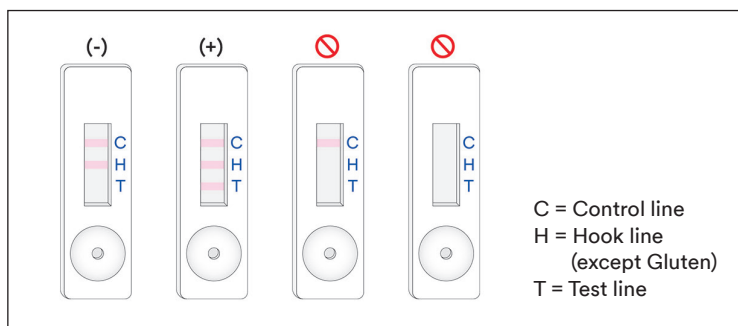Any reading after 12 minutes from the initial application of the sample into the Neogen Soy Protein LFD should be considered invalid. A reading at this time cannot be interpreted and can lead to erroneous results.

## Troubleshooting

1. Sample fails to migrate across the strip within the first 5 minutes after application of the sample in the Neogen Soy Protein LFD.  
The sample may be too viscous and needs to be centrifuged if this was not already done during the preparation of the sample. If sample was already centrifuged, then preparing a 1:1 dilution with the Neogen Extraction Buffer may be necessary. (Note: This may reduce the sensitivity to ~4 ppm for some matrices.)
2. A red dot appears on the test line but remainder of test line does not change color.  
Sample particulate may have passed around the filter in the cassette, simply re-run the sample by taking a new Neogen Soy Protein LFD from the kit and repeat the test.
3. The pH of extracted sample should be between 5 and 10. If pH is outside this range, further dilution may be required (i.e., prepare a 1:1 dilution with 100  $\mu$ L of the extracted sample and 100  $\mu$ L of Neogen Extraction Buffer. This may reduce the sensitivity to ~4 ppm for some matrices.)

If you have questions about specific applications or procedures, please contact your Neogen Food Safety representative or distributor.

## Minimum Performance Characteristics

Lowest Limit of Detection <sup>(a)</sup>	2 ppm
Upper Limit of Detection	~2 mg/mL

- <sup>(a)</sup> The lowest limit of detection is defined as the lowest concentration of the allergen in a test sample that can be distinguished from a true blank sample at a specified probability level<sup>1</sup>.

## References

1. Abbott, M., Hayward, S., Ross, W., Godefroy, S.B., Ulberth, F., Van Hengel, A. J., Roberts, J., Akiyama, H., Popping, B., Yeung, J.M., Wehling, P., Taylor, S., Poms, R.E., and Delahaut, P. (2010). Appendix M: Validation Procedures for Quantitative Food Allergen ELISA Methods: Community Guidance and Best Practices. *J. AOAC Int.* 93, 442-450.

## Explanation of Symbols

[info.neogen.com/symbols](http://info.neogen.com/symbols)

# Neogen Food Safety

## Neogen Corporation

620 Leshar Place  
Lansing, MI 48912 USA  
Neogen.com

## Neogen Europe Ltd.

The Dairy School  
Auchincruive  
Ayr, KA6 5HU  
Scotland, UK

## Neogen Ireland, Ltd.

Bray Business Park, Bray  
Co. Wicklow  
A98YV29, Ireland

## Neogen Corporation

620 Leshar Place Lansing, MI 48912 USA  
www.neogen.com

Neogen is a trademark of Neogen Corporation.  
© Neogen Corporation 2024.  
Unauthorized use prohibited. All rights reserved.

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